



DURAN® BOTTLE CAPS & CONNECTION SYSTEMS

THE **PERFECT** MATCH – FROM STANDARD TO SPECIAL

In scientific procedures, nothing should be left to chance. Factors such as functionality, efficiency and cost-effectiveness in laboratories can be affected by apparently simple elements such as bottle caps and connection systems.

DURAN® offers a comprehensive product range of precision closure systems, starting with high-quality standard caps suitable for numerous scientific applications.

Our range also includes identification aids, special closures for sterilization processes, caps made of high performance materials for demanding applications, connection systems, security caps and transportation accessories. DURAN® caps and connection systems fit most DIN standard laboratory bottles made by other manufacturers.





Find out more:

www.duran-bottle-system.com



THE EVOLUTION OF THE **DURAN®** BOTTLE SYSTEM

Starting with the ORIGINAL in 1972, the DURAN® laboratory bottles have continued to evolve to meet the diverse needs and applications of scientists and researchers around the world. The borosilicate glass 3.3 bottles are now available in the most comprehensive range of volumes and formats.



FOR ALL YOUR NEEDS

DURAN® BOTTLE ACCESSORIES

A range of dedicated accessories that work perfectly with the DURAN® laboratory bottles. DURAN® accessories enable practical color coding, product labelling and identification. Most are reusable and long lasting just like the DURAN® bottles.





Click and Carry! Our latest innovation – a modular bottle carrying system that focuses on the safety and convenience of users in the laboratory.

Find out more on page 10

SELECT THE BEST COMPONENTS FOR YOUR NEEDS

STANDARD



These reliable components are made of foodgrade polypropylene (PP) and are suitable for many laboratory applications. Our screw caps and pouring rings are fully compatible with DIN GL bottle threads.

HIGH PERFORMANCE



For more demanding processing conditions, or when the product requires a higher level of protection from the environment or contamination. Ideal for demanding applications in the pharmaceutical industry.

STEAM STERILIZATION



Membrane caps are recommended when sterilizing media. These components improve the safety and ease-of-use of the laboratory bottles during autoclaving or hot-air sterilization, and protect the contents from contamination.

CONNECTION



The robust and durable DURAN® tubing connection systems allow the easy transfer of liquids under controlled conditions, reducing the risk of contamination. Connector elements are also available for probes and stirrers.

SECURITY AND TRANSPORTATION



When extra protection is required to ensure the safety of the contents and users during storage or transport. For example, by blocking ultraviolet light, ensuring contents aren't tampered with or making transport easier and safer.

IDENTIFICATION



The accurate identification of samples and reagents is integral to the scientific process. These components provide simple color coding or labelling of laboratory bottles.

NECK SIZE

GL25

NECK SIZE

GL32

Original GL Screw Cap Blue, PP, Linerless

ANSI/SLAS Microplate Holder for 10 or 25 mL bottles



High Temperature Screw Cap Red, PBT, Cap liner



Aperture Screw Cap Red, PBT



Premium Cap

PURE Premium Screw Cap Colorless, PFA, Cap liner



Stopper Grey, Bromobutyl rubber



Original
GL Screw Cap
Blue, PP, Linerless



Original Pouring Ring Blue, PP



High Temperature Screw Cap Red, PBT, Cap liner



Pouring Ring Aperture Red, ETFE Screw Ca





Screw Cap Membrane Vented Blue, PP



Multiport Connector Cap 2-, 3- or 4-ports, PTFE with stainless steel connectors



Torque wrench socket

fitting for PP Cap*



Torque wrench socket fitting for PBT Cap*



Torque wrench socket fitting for PREMIUM Cap*



Screw Cap Membrane Vented Blue, PP



Multiport Connector Cap 2-, 3- or 4-ports, PTFE with stainless steel connectors



Torque wrench socket fitting for PP Cap*



Torque wrench socket fitting for PBT Cap*



GL 45

STANDARD

PERFORMANCE







Original GL Screw Cap Blue, PP, Linerless



Original Pouring Ring Blue, PP



YOUTILITY® Screw Cap Cyan, PP, Linerless



YOUTILITY® Pouring Ring Cyan, PP



High Temperature Screw Cap • Red, PBT, Cap liner



Pouring Ring Red, ETFE



YOUTILITY® Screw Cap **Aperture** USP Class VI compliant **Screw Cap** Red, PBT White, PP, Linerless



Premium Cap Colorless, PFA, Cap liner



Premium Pouring Ring Colorless, PFA



Screw Cap Membrane Vented Blue, PP



YOUTILITY® Screw Cap Membrane Vented Cyan, PP



Multiport Connector Cap 1-port, Stainless Steel



Multiport Connector Cap 2- or 3-ports,



Multiport Connector Cap 2-, 3-, 4-or 5-ports, PTFE, Stainless Steel Connectors for narrow Tubing



Multiport Connector Cap 4-ports, Stainless Steel



Screw Cap 2 Hose Connectors Blue, PP



Multiport **Connector Cap** 2- or 3-ports, Blue, PP



Multiport Connector Cap **HPLC 3-port** Blue, PTFE



Tamper Evident Screw Caps Yellow or Red, PP, Linerless or Cap liner



Torque wrench socket fitting for PP Cap*



Torque wrench socket fitting for PBT Cap*



Torque wrench socket fitting for PREMIUM Cap*



Original GL Screw Cap Yellow, Green or Grey,



Pouring Ring Yellow, Green or Grey, PP



YOUTILITY® Screw Cap. Transparent, PP, Linerless



YOUTILITY® Bottle Tags Various colors, Silicone





By using an innovative 3D metal printing process, we have optimized the design of the Multiport Connector Cap for liquid transfer applications.







Aluminium

Aluminium Screw Cap Open top



Over-cap Black, Phenolic



Stopper Grey, Bromobutyl rubber



Septum Colorless, PTFE coated Silicone



Multiport **Connector Cap HPLC 4-port** Blue, PP



Screw Cap for pH Sensor Red, PBT



Replacement Cap Liner PTFE coated Silicone



Septum Colorless, Silicone



Stirred Reactor



Screw Cap with Temperature Probe



Screw Cap with Temperature Probe





YOUTILITY® Printable Labels, PS






The DURAN® bromobutyl rubber stopper provides a gas tight seal for all GL 45 laboratory bottles. Bromobutyl rubber is essentially impermeable to most gases, and provides a controlled environment inside the glass bottle for oxygen sensitive materials.

NECK SIZE

GL 56

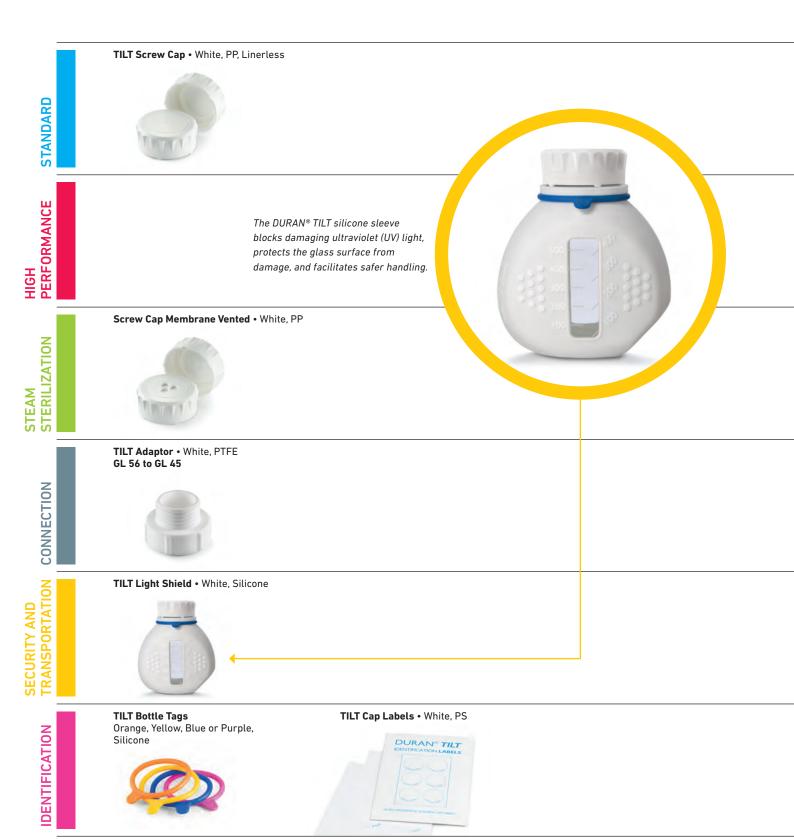

NECK SIZE

GLS 80[®]

In addition to the benefits of a larger mouth, such as hassle-free filling, emptying, and cleaning, the GLS 80® caps have a special fast closing thread.

Screw Cap • Blue, PP, Linerless



Pouring Ring • Blue, PP



Screw Cap • White, PSU, Cap liner



PURE Screw Cap • White, PSU, Cap liner



Pouring Ring • White, PTFE



Screw Cap Membrane Vented • Blue, PP



Multiport Connection Cap • 4-ports, Blue, PP



 $\textbf{Stirred Reactor} \bullet \text{Blue and gray, PP}$



Connection Cap System for Overhead Mixer • PP





The DURAN® Stirred Reactor is ideal for the widest possible range of mixing processes in laboratories, e.g. mixing liquids or dissolving solids. The variable stirring shaft can be used for 500 ml, 1000 ml and 2000 ml DURAN® GSL 80® laboratory bottles and can be driven with standard magnetic stirrers.

DURAN® GL 45 BOTTLE CARRYING SYSTEM

The modular Bottle Carrying System fits standard laboratory bottles that have a GL 45 neck finish and is compatible with four bottle sizes: 2, 5, 10 and 20 litres. It is ergonomically designed, and offers health & safety benefits to the user and allows easy and safe transport of bottles around the laboratory or production facility. It makes pouring liquids from the larger bottles more comfortable and safer. All materials are fully autoclavable.

Carrying Handle

- Carrying handle is suitable for bare or gloved hands
- Handle strap is flexible enough to allow pouring

Tightening Strap

- With indicator arrow
- End modified to prevent unthreading

Pouring Handle

- Ergonomically shaped
- Anti-slip grip

Silicone Support Base

- Structural support for base of larger 5, 10 and 20 litre bottles
- Protection of glass bottle base from impact during lifting or handling
- Acts as a spacer to protect glass from contact with Metal Dolly
- Moulded-in holes minimize liquid traps and allow maximum steam penetration during autoclaving



GL 45 Bottle Collar

- The collar is securely attached by pushing down over the GL 45 thread with 4 teeth that lock under transfer ring.
- The collar is only removable by cutting so attachment is 100% secure in use
- A cut line indicator for removal of collar when required

Horizontal Straps

- Durable stitched construction
- Made from "anti-slip" material



Lot Specific Retrace Code for traceability and certification



Optional Extra: Metal Dolly

- Stainless steel dolly for 10 or 20 litre bottles (2 sizes only)
- Fully compatible with carrying system
- Easy to clean, crevice-free construction
- Optional handle is available for easier handling



FOR EASY AND SAFE CARRYING & POURING



Easy and safe carrying using handle

• Handle attaches at the neck and fully supports the base



Safe and controlled pouring using handle

 Safer pouring of liquids using a second lower "Pouring Handle" located on the side strap (Available on the 5 litre carriers only)

Fast and secure tightening

of the polypropylene straps with two side strap adjusters for the larger 5, 10 and 20 litre bottles



2 L bottle carrier in 4 colors


for color coding









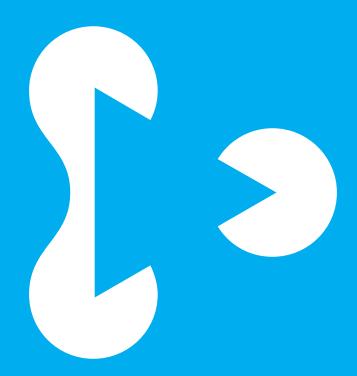

Excellence in your hands



Phone: +49 6131 - 1445 4131 Fax: +49 6131 - 1445 4016 sales@DWK.com www.DWK.com

